

## Updating or Installing Insync

InSync can be downloaded from the following location:


<https://www.escription-one.com.au/Software/InSync>

1. Click on the download hyperlink for version 6.11 as noted in the below image. Your download should start.


**InSync Software Download**

**!** Support is only available for products running on Windows 10 and later operating systems.


One of the steps of installing and configuring your recorder will direct you to download the InSync software, found below.

Version	Install Link	Manual	Version History
6	<a href="#">6.11.20329</a>		<a href="#">Version History</a>


2. When downloaded, **Run** the file.
3. You will be presented with the InSync setup wizard. Click **Next** to continue with the installation.
4. Read the license agreement, then select the '**I accept the terms in the license agreement**' checkbox. Click the **Next** button.
5. Unless otherwise specified by your IT department, choose the option to **Install for all users of this machine**. Then click the **Next** button.

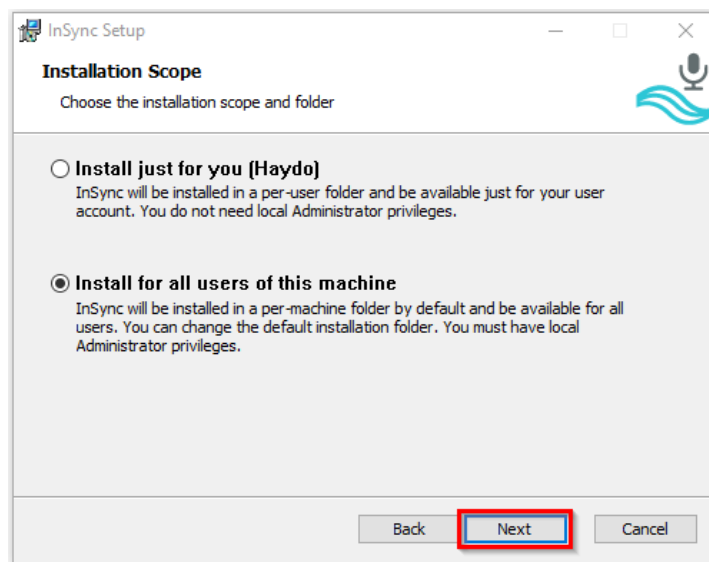


6. Unless otherwise specified by your IT department, do not change the default installation path. Click the **Next** button.

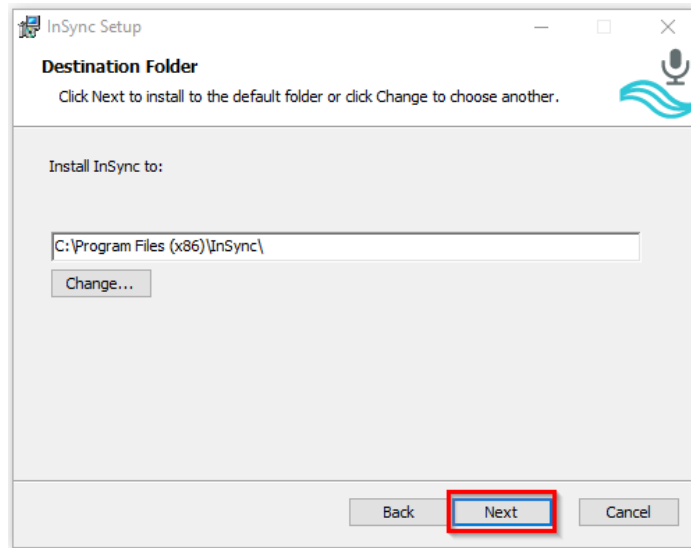


7. Click the **Install** button. InSync will then install to your system. You may need to accept a UAC prompt to allow the install.

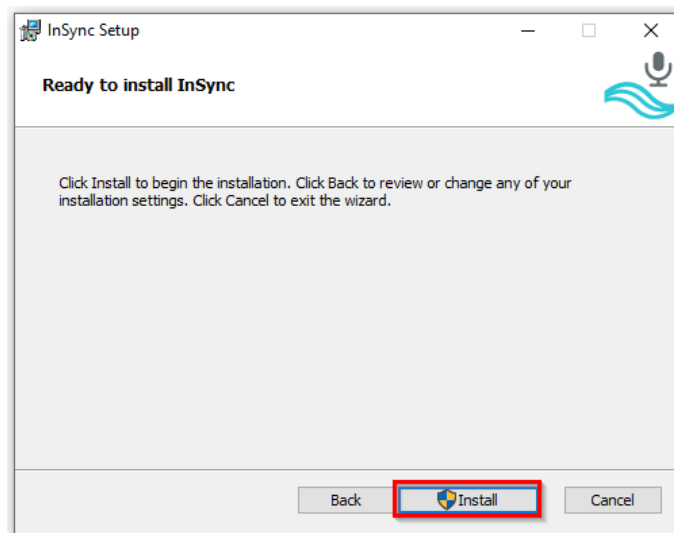

8. Select check boxes in the Finish window if desired. **'Start InSync with Windows'** will start the application with Windows. Please check this box if you would like InSync to start automatically with windows. **Click Finish** to complete the installation.

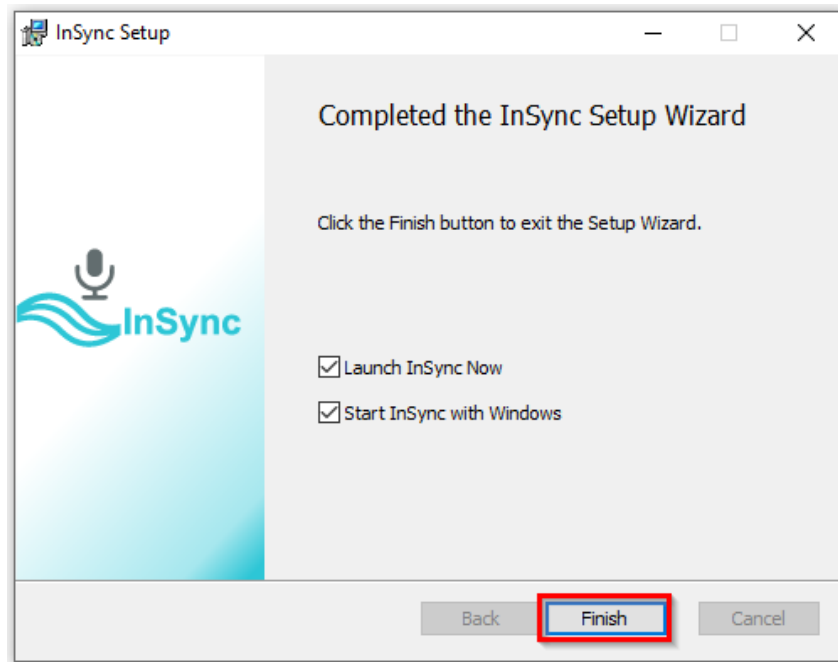

9. When InSync is launched please click on Change Region and select Australia as the region.

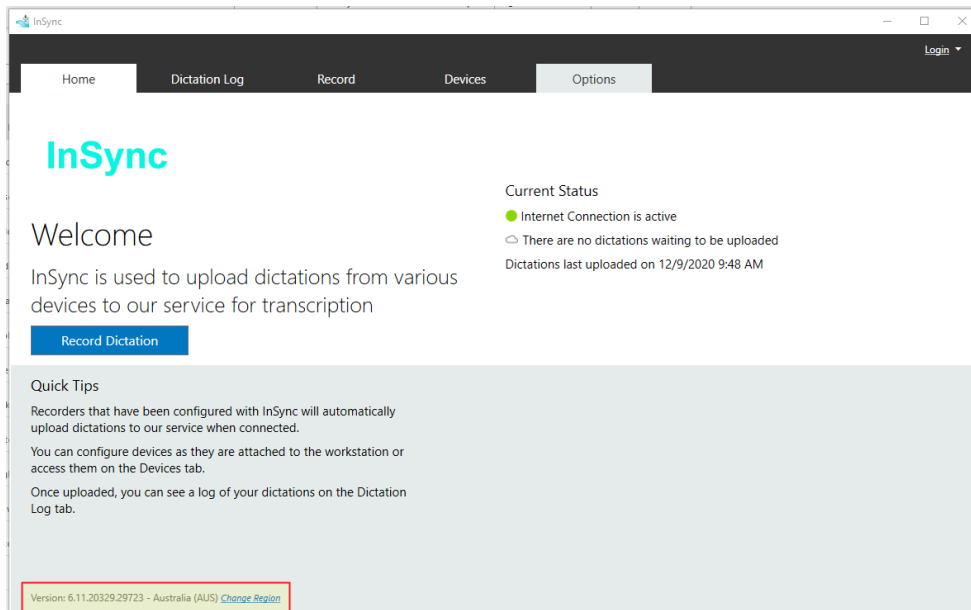

Click Yes on the "InSync must restart ..... " message to complete the upgrade.

